Metropolitan Museum of Art Gas Chromatography- Mass Spectrometry (GC-MS) Results from Material Analysis

This document includes (1) a mass spectrum and (2) the volatile organic compounds (VOCs) emitted from samples using GC-MS analysis. The data is not interpreted; however, several classes of chemicals are highlighted because they are potential risks for artwork in an enclosed environment. A basic key, provided below, indicates those classes. The amount of each chemical identified has not been determined; similarly, it is not known how much of each chemical is necessary to do damage to art. Finally, peaks may be present that are the result of the sample adsorbing chemicals from the air and reemitting them during testing rather than being inherent to the sample. Research is ongoing to determine specifically which chemicals and amounts are required to negatively affect artifacts.

Highlighted data:

Pink – chemicals currently known to be hazardous to art

Green – amines; can raise the pH, are suspected to react with acids and may form crystals in an enclosed environment

Yellow – chemicals of the following type, which may be hazardous to art:

Acids – lower the pH, corrosive to metals, degrade organic materials

Aldehydes – can convert to acids with heat or exposure to UV light

Esters – can hydrolyze into acids with heat and humidity

Sulfur-containing compounds – known to tarnish and corrode some metals

Halogenated compounds – can become reactive with exposure to heat and UV light

Nitrogen-containing, not amine – can react with other off-gassed chemicals


Alkynes – can become reactive when exposed to heat or UV light

Sample: Creation Baumann Velling III 0366 cotton velvet fabric in powder blue

Oddy test result: Temporary Date collected: 8/18/2018

Technique used: SPME with a PDMS/DVB fiber; Agilent 7890B GC and 5977B MS fitted with a GL Sciences OPTIC-4 multimode inlet and LEAP PAL RTC autosampler; Pre-heated at 60°C for 20 minutes; fiber exposure at 60°C for 20 minutes; sample injected into 220°C inlet and cryo-trapped for 2 min at -15°C; GC ramped from 35°C to 250 °C at 10°C/min. Data analyzed in Masshunter Qualitative Analysis. Deconvoluted data with > 85% match with a NIST 17.0 or Wiley 9 library are reported.

VOCs not highlighted are because they were also observed in blanks: (1) 11.6 min: 2-methyl-, 2,2-dimethyl-1-(2-hydroxyl-1-methylethyl) propyl ester propanoic acid; (3) 11.9 min: 2-methyl-, 3-hydroxyl-2,4,4-trimethylpentyl ester propanoic acid

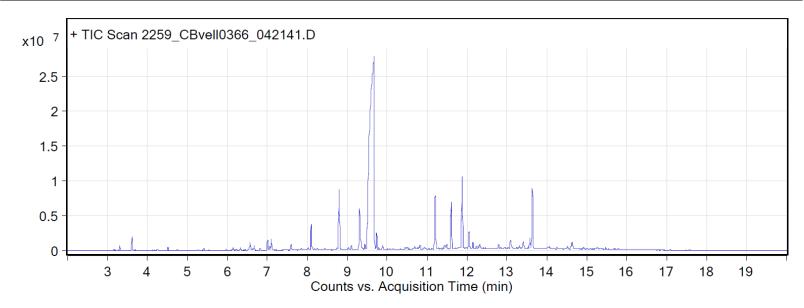

Compound Table

RT	Score (Lib)	Area	Name	Formula
6.56	96.23	1468172	Cyclotetrasiloxane, octamethyl-	C8H24O4Si4
7.01	93.85	2042529	1-Hexanol, 2-ethyl-	C8H18O
7.1	95.16	2299092	Benzyl Alcohol	C7H8O
8.09	97.44	5474359	Nonanal	C9H18O
8.8	94.87	13359851	Cyclopentasiloxane, decamethyl-	C10H30O5Si5
9.31	95.12	9400797	Ethanol, 2-(2-butoxyethoxy)-	C8H18O3
9.73	97.32	2749012	1,3-Dimethyl-2,4,5-trioxoimidazolidine	C5H6N2O3
11.2	96.01	12011795	Cyclohexasiloxane, dodecamethyl-	C12H36O6Si6
11.6	89.97	10445037	Propanoic acid, 2-methyl-, 2,2-dimethyl- 1-(2-hydroxy-1-methylethyl)propyl ester	C12H24O3
11.87	93.74	16093824	Propanoic acid, 2-methyl-, 3-hydroxy- 2,4,4-trimethylpentyl ester	C12H24O3
12.05	96.18	4221907	1-Tetradecene	C14H28
12.15	94.5	1394176	Tetradecane	C14H30
13.57	92.13	1819933	Phenol, 2,4-bis(1,1-dimethylethyl)-	C14H22O
13.64	96.75	11/85408	Phenol, 2,6-bis(1,1-dimethylethyl)-4- methyl-	C15H24O
14.63	93.1	1956500	PENTAN-1,3-DIOLDIISOBUTYRATE, 2,2,4-TRIMETHYL-	C16H30O4